

Japan Industry-Government- Excellence Award 2019 Award winner's product

The keyword is "20 seconds".

"The Kitsutsuki Motion Contra" for root canal enlargement of startle debuts as the top since fall of 2018.

Until now, the contra canal for enlarging the canal enlarges the canal by low-speed rotational movement. However, root canal enlargement due to slow rotation is time consuming and has the disadvantage that the file is prone to fracture due to "twisting" caused by rotation. The root canal enlargement contractor "The Kitsutsuki Motion Contra" we introduced this time has made it possible to expand the root canal in an extremely speedy and safe manner.

Please try a world of surprise as never before.

● Vertical stroke not rotation

Instead of rotating the attached file as in the past, the "rooting motion" that moves the file up and down like a woodpecker efficiently cuts and enlarges the root canal. Also, you can use the file as it is.

● Completed in 20 seconds

The biggest feature of this contra is to expand the canal in an extremely speedy and very easy way. In the clinical trial, it took about 20 seconds, and the extraction was completed in about 9 seconds.

● Relief of the patient's pain

Since it ends in a short time vertical enlargement, lateral enlargement, and wall grinding, the patient's pain is greatly reduced.

● File corruption is drastically reduced

As mentioned earlier, files that are elongated needle-like have the disadvantage of causing metal fatigue due to "twisting" by rotation and breaking, but this woodpecker motion controller is a structure that moves files up and down. Even in the case of a bent root canal, breakage is unlikely to occur, and since there is no eccentric enlargement, enlargement work can be performed safely.

● High-performance, high-quality handmade, so because Made in Japan

The woodpecker motion contras are a collection of fine manufacturing technologies that Japan boasts. And because the surface is given a high-grade plating coating, it is resistant to scratches and has become a product with excellent durability. The 2-year product warranty is an indication of the excellence of Made in Japan



Authorized Distributor (Overseas): Miletusthales Co., Ltd /MS. R. HIDAHA

WWW.MILETUSTHALES.ORG

Classification: Contra-angle for root canal enlargement

(Specific maintenance management medical equipment / approval number
230ALBZX00028000)

Specifications

Head: Standard head

Light: no light

Water injection function: None

Gear ratio: deceleration 4: 1

Head movement: vertical stroke movement + free rotation

Mounting file: H file (Manny super file recommended)

Head size: Ø11 x H18.5 mm (including the rear screw)

Number of revolutions: 20,000 min-1


(recommended number of revolutions 8,000 rpm or less)

Compatible micromotor: can be mounted to an ISO standard micromotor

Sterilization: autoclavable (135 ° C)

Warranty period: 2 years from date of purchase

made in Japan

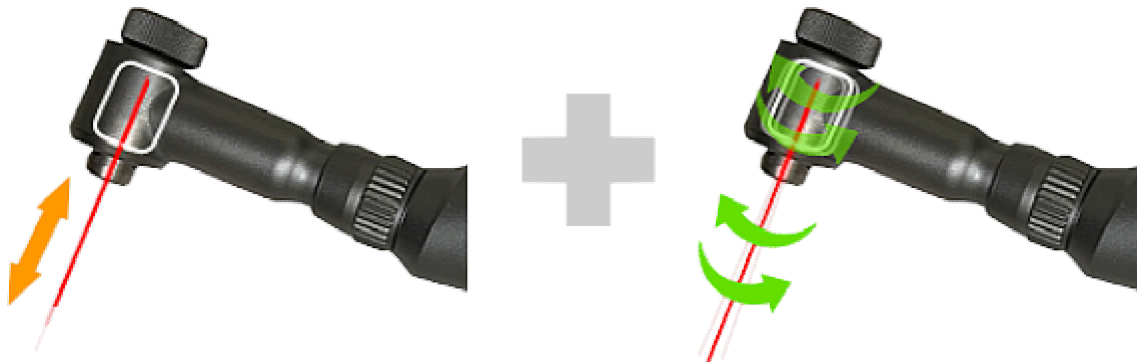

Authorized Distributor (Overseas): Miletusthales Co., Ltd /MS. R. HIDAKA

WWW.MILETUSTHALES.ORG